

Date: October 15, 2009


TN MEK101509

Subject: Replacement of MEK Rub Test Machine Safety Cover

Models: All

Scope: Procedure for installation of replacement of safety cover

1. Remove the old cover, keep mounting hardware and discard cover.
2. Remove protective plastic film from the new cover.

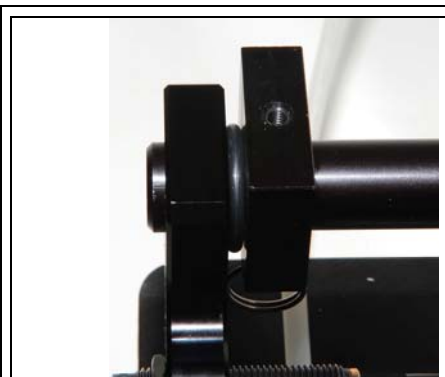

3. Align mounting blocks and ensure that they are pushed out towards the cylinder mounting brackets. This will eliminate any side play in the cover.



4. Center the cover on the machine and align the front edge of the new cover to the front edge of the base plate as shown.



5. Mark the position of the mounting holes, remove the cover and drill the marked positions 5mm or 13/64" \varnothing .

To prevent binding or cracking when drilling the holes place a piece of wood under the Lexan and drill through the cover into the wood.



6. Install the mounting hardware as shown